

USER'S
MANUAL

FLPG Plus




CONTENTS

Before the first use.....	3
Product Overview.....	4
Quick Start Guide.....	5
Output Capacity.....	8
Specifications.....	9
Technical Support.....	10
Notes.....	11

BEFORE THE FIRST USE

Remove the transport screws under the FLPG Plus.

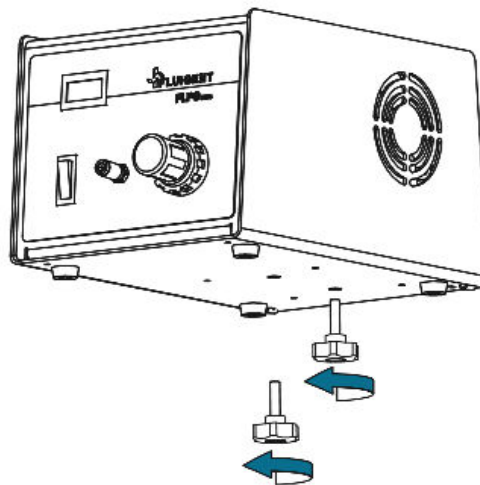


These screws are used to protect the compressor inside during the transport.


Please keep these screws in case of future transport.

PRODUCT OVERVIEW


The pressure value displayed is in units of bar (1 bar = 14.5 psi)

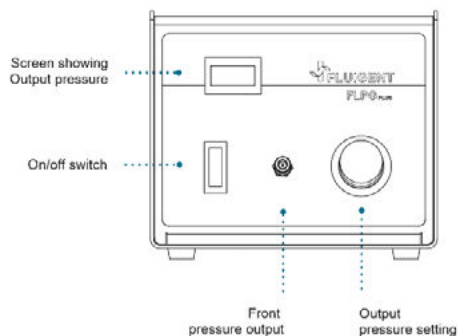

The rear pressure output delivers the same pressure as the front output. One or the other can be used for better tubing management, or both at the same time to pressurize two pressure controllers.

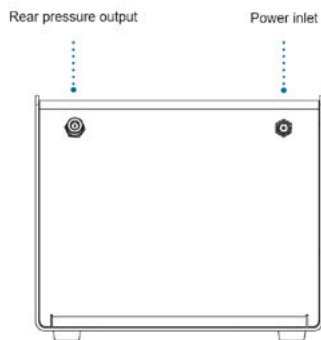


QUICK START GUIDE

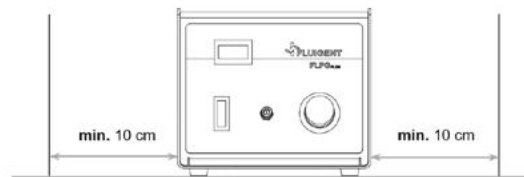

Please, always place the FLPG Plus on a horizontal and stable surface. Keep clear for at least 10 cm at sides of FLPG Plus. This allows for ventilation and the air drying system to function properly.

START USING THE FLPG PLUS

- 1 Supply and power ON the FLPG Plus
- 2 Pull out the pressure regulator knob to unlock it

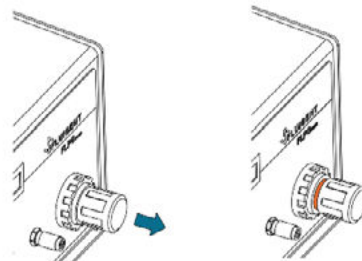

Locked

Unlocked (orange ring is visible)

- 3 Turn the knob to change the output pressure (from 0.1 bar to 2.30 bar) in order to meet the supply requirement of the pressure controller channel range.

Different channel ranges require different pressure supply. The following table shows the cases of Fluigent's pressure controllers.


Flow EZ™ (channel range)	2000 mbar	1000 mbar	345 mbar	69 mbar
Required pressure supply	2.1 bar	1.1 bar	1.1 bar	0.15 bar
MFCST™ series (channel range)	2000 mbar	1000 mbar	345 mbar	69 mbar
Required pressure supply	2.3* bar	1.3 bar	0.8 bar	0.5 bar

*The MCFS 2000 mbar range channel is compatible with the 2.3 bar supply pressure, despite its nominal supply pressure of 2.4 bar


The FLPG Plus is designed to work at 2.30 bar or below.

Do not operate above this value, or the output pressure can become unstable. This may damage the FLPG Plus and the connected devices.

- 4 Connect the pressure controller's pressure inlet to the FLPG Plus' pressure output (front or/and rear)

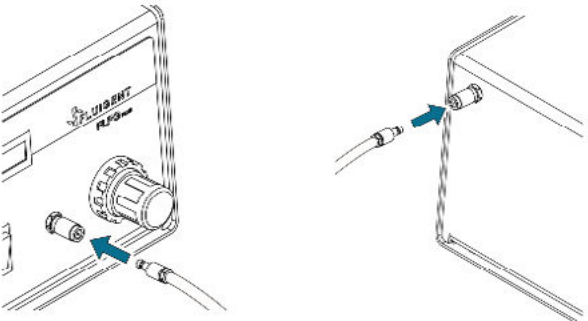

- 5 Power OFF the FLPG Plus after use. Doing so will extend service life of the compressor

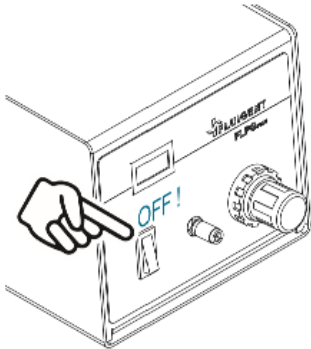

OUTPUT CAPACITY

The FLPG Plus' output capacity is the number of channels it can drive simultaneously. This capacity depends on the output pressure set with the pressure setting regulator.

When set to 2.1 bar to drive the Flow EZ™ 2000 mbar channels (or 2.3 bar to drive MFCS 2000 mbar channels), the FLPG Plus can drive up to 4 channels at the same time.

When set to 1.1 bar to drive the Flow EZ™ 1000 mbar channels (or 1.3 bar to drive MFCS 1000 mbar channels), the FLPG Plus can drive up to 8 channels at the same time.

When set to lower pressures to drive lower range channels (see tables on page 6 for more information), the FLPG Plus can drive up to 8 channels at the same time.

Do not connect more channels than the numbers shown above. Otherwise, the output pressure may drop and become unstable, leading to possible damage of the FLPG Plus and the connected pressure controller.

SPECIFICATIONS

Product Number: FLPG003

Output pressure: From 0.1 bar to 2.3 bar
2 bar: up to 8 channels
1 bar: up to 16 channels

Voltage: 24V DC

Current: 1.6 A

Dimensions: 179 x 154 x 264 mm (W x H x D)

Weight: 5.9 kg

Environment: Indoor clean environment
15 -25 °C
RH < 65% (no condensation)

Lifetime: 10 000h**

* Higher output capacity version exist, consult Fluigent sales representative

**Maintenance is recommended after this duration

Page 10 of 10

Still have questions? E-mail us at:
support@fluigent.com

Still have questions? E-mail us at:
support@fluigent.com

Or call our technical support team directly



Fluigent S.A.S +33 1 77 01 82 65

For a fully detailed FAQ for all Fluigent products, please visit:



<http://www.fluigent.com/faqs/>

Interested in Fluigent products ?

To view the [complete Fluigent product line](#) and [application notes](#):



<http://www.fluigent.com>

For commercial requests, please e-mail:




contact@fluigent.com or your local office

For [tutorial videos](#) about the [LineUp™ series](#), please visit [Fluigent on YouTube](#)



Fluigent

Page 10 of 10

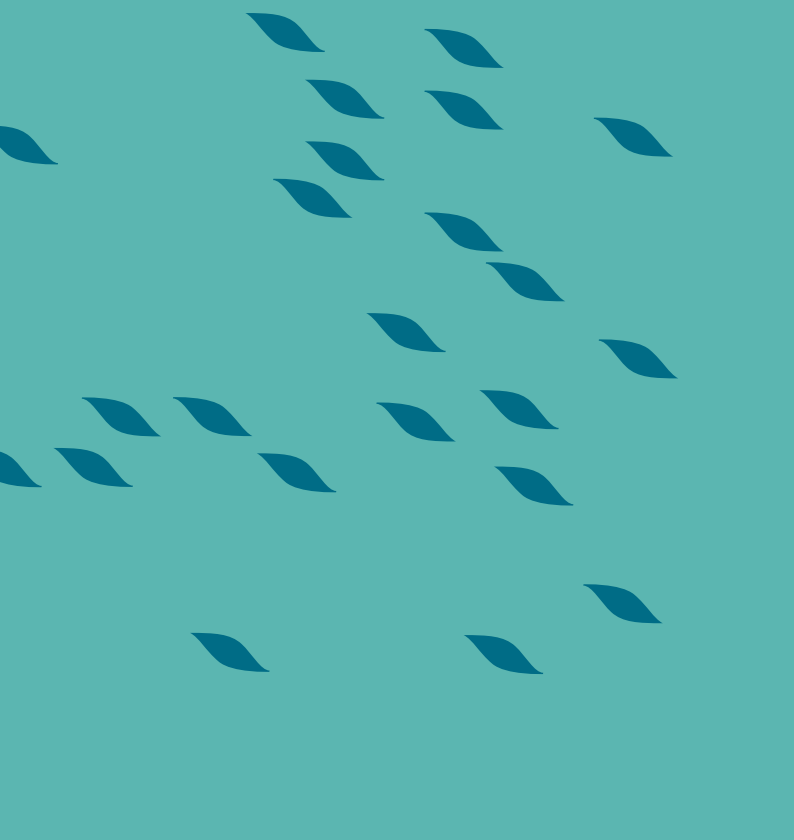

**VERSION
2026**

